


SOLUTIONS FOR YOUR LABORATORY


VERDER SCIENTIFIC IoT APP QUICKSTART GUIDE

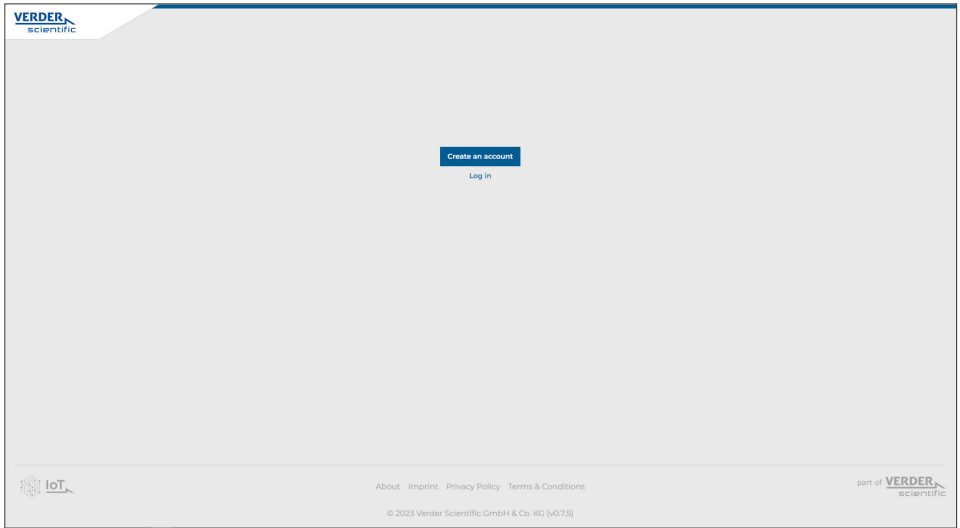

1. GO TO [IOT.VERDER-SCIENTIFIC.COM](https://www.iot.verder-scientific.com), CLICK THE „CREATE AN ACCOUNT“ BUTTON

This is the first step in creating an account on the platform.

You will need to locate and click on the button that says „create an account“.

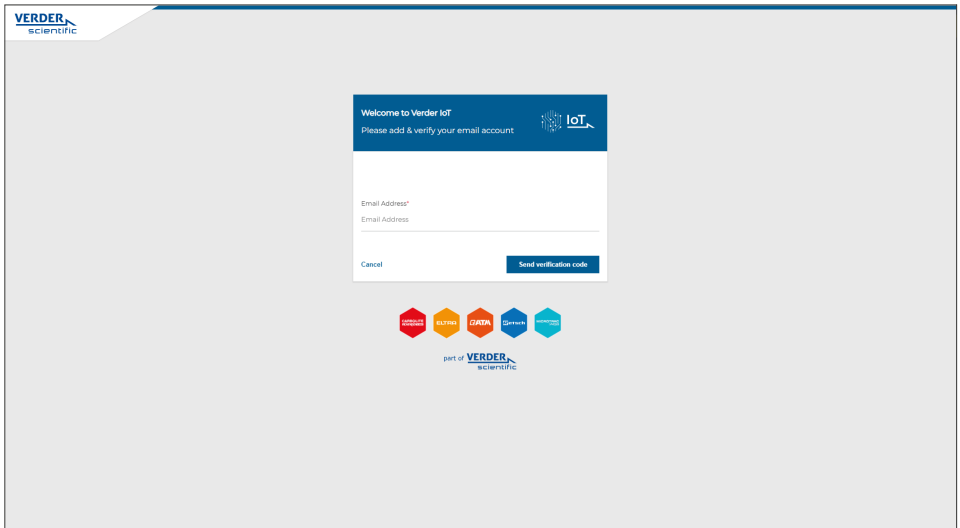

2. INSERT YOUR EMAIL

After clicking the create an account button, you will be prompted to enter your email address.

This is the email address that will be associated with your account on the platform.

Welcome to Verder IoT
Please add & verify your email account

Verification code has been sent to your inbox. Please copy it to the input box below.

Email Address*
e@p.lenc@verder.com

Verification Code*
154003

Cancel Verify code Send new code

part of **VERDER** scientific

3. VERIFY YOUR EMAIL BY INSERTING THE CODE SENT TO YOUR EMAIL

Once you have entered your email address, a verification code will be sent to that email (please check your spam folder). You will need to retrieve the code from your email and enter it into the platform to verify your email address.

Welcome to Verder IoT
Please fill in your account information

New Password*
Show Password

Confirm New Password*
Confirm New Password

I am a Customer*
Customer

Personal Information

Name*
Surname

Email Address*
Email Address

Business Information

Company*
Address

Account Information


Username*
Password*
Email Address*

Create

part of **VERDER** scientific


4. SET A PASSWORD, CREATE A TENANT AND A LAB, AND FILL IN YOUR BUSINESS PROFILE INFORMATION

After verifying your email, you will be prompted to set a password for your account. You will also need to create a tenant and a lab, and fill in your business profile information. This information will be used to set up your account on the platform.

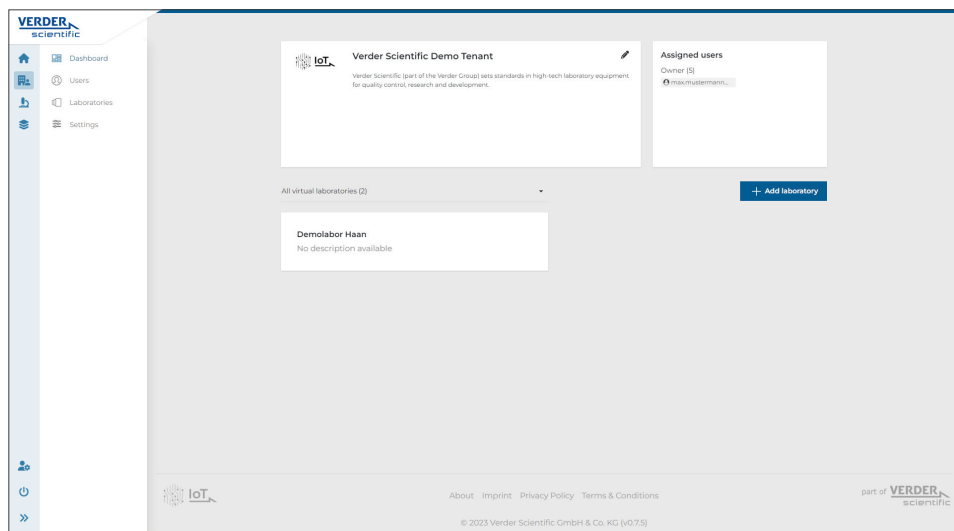

5. NAVIGATE TO THE LAB BY CLICKING ON THE CARD WITH YOUR LABS NAME


Once your account is set up, you can navigate to your lab by clicking on the card with your labs name. This will take you to the lab where you can begin using the platform

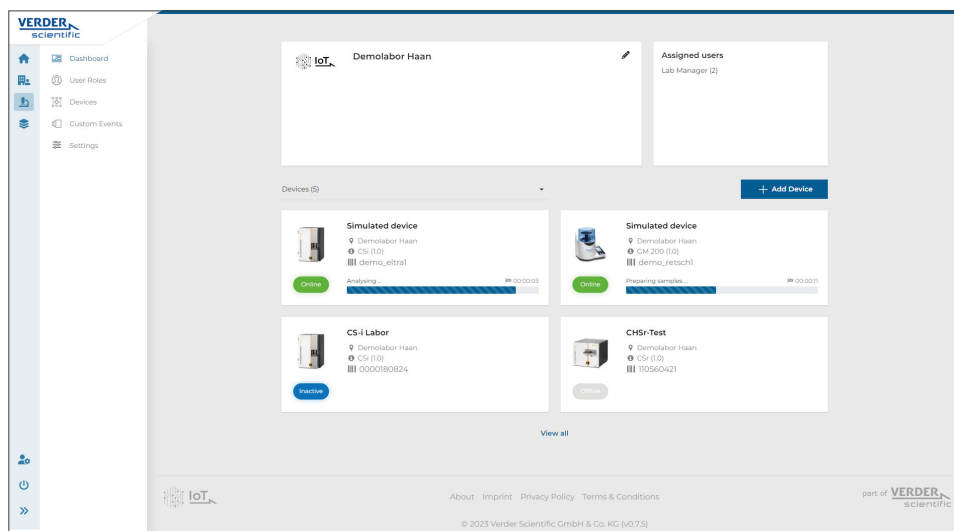

6. CLICK THE ADD DEVICE BUTTON


In the lab, you will see an „add device“ button.

Click on this button to begin the process of adding a device to the platform

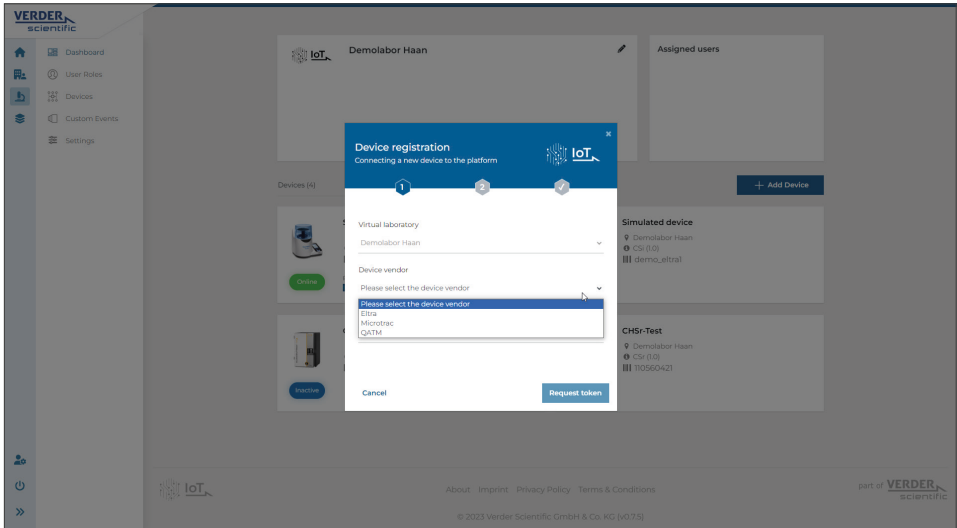


7. SELECT THE VENDOR OF YOUR DEVICE

After clicking the add device button, you will be prompted to select the vendor of your device. This is the company that manufactured your device.

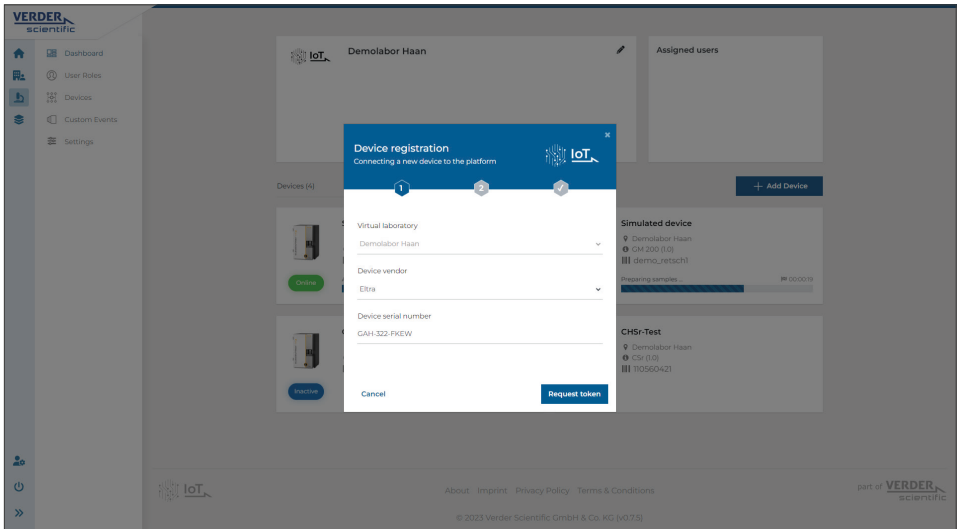


8. ADD THE DEVICE SERIAL NUMBER

Once you have selected the vendor of your device, you will need to enter the device's serial number. This is a unique identifier for your device.

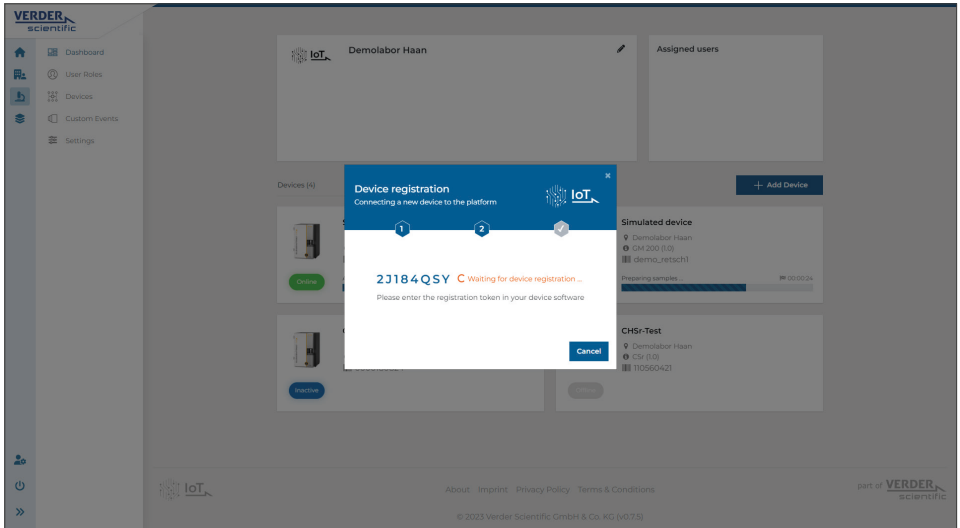

9. Enter the code displayed on the platform into your device software

After entering the device serial number, a code will be displayed on the platform. You will need to enter this code into your device's software to complete the process of adding the device to the platform.

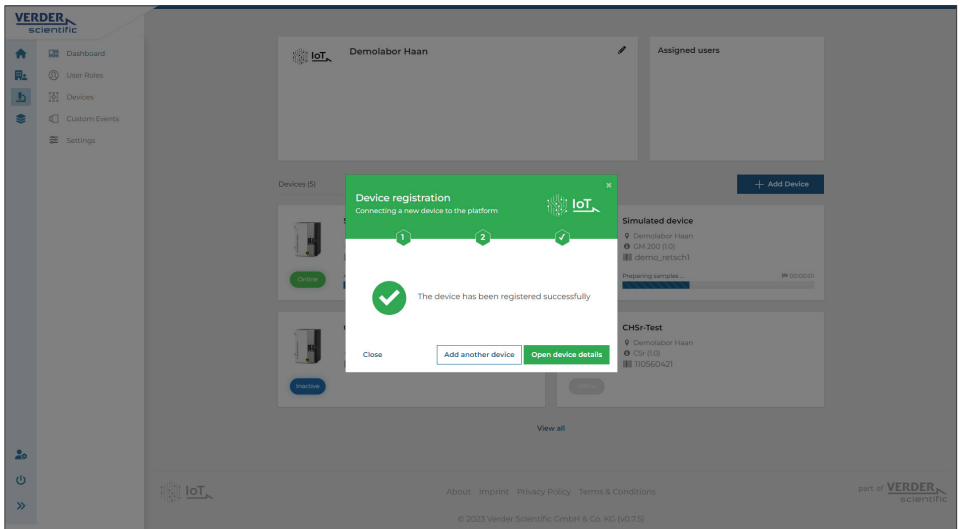

10. The operation has been successful, and you are now onboarded on the platform

If all of the previous steps have been completed successfully, you will receive a message indicating that the operation was successful and that you are now onboarded on the platform. You can now begin using the platform with your device.

NOTES



VERDER SCIENTIFIC GmbH & Co. KG

Retsch-Allee 1-5

42781 Haan · Germany

info@verder-scientific.com

www.verder-scientific.com